FeverVue™ Digital Thermometer


Item 15-905-000

Thank you for purchasing a HealthSmart Digital Thermometer. Please read the following instructions to achieve the most accurate temperatures and safe operation. With proper care and use, your thermometer will provide you with reliable readings. It is enhanced with a beeper, memory and fever alarm. The following instructions will cover your thermometer's features in further detail.

What is a "Normal" Temperature?

Although the generally accepted "normal" temperature reading is 98.6°F(37°), temperature readings can vary from 97.0°F(36.1°C) to 99.0°F(37.2°C) and still may be considered "normal". Variations in temperature can be attributed to activities such as exercise, smoking, eating and drinking. Even time of day may influence your temperature. For example, your temperature is lower in the morning than in the afternoon. Other variations may be due to the location of the temperature reading. While the oral temperature follows the guidelines previously described, a rectal temperature is generally 1.0°F higher. Conversely, an axillary (under the arm) temperature will be 1.0°F lower.

Features of Your Digital Thermometer

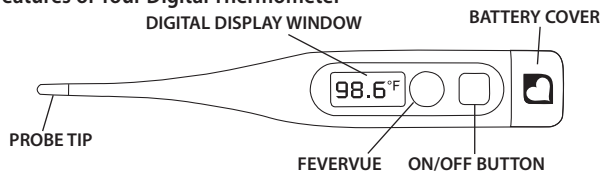Cleaning Instructions

Before and after each use, clean the thermometer with a soft cloth and 70% isopropyl alcohol diluted with water, or wash with COLD soapy water. Wipe again with a clean, water-dampened cloth and let dry.

- The unit must never be sterilized by boiling, gas or steam autoclaves.
- The use of any other solution or cleaning method beyond the above mentioned may result in reliability problems with the thermometer and void your warranty.

How to Take Your Temperature

NOTE: Measurement time may vary based on the selected site.

- Clean and disinfect the thermometer according to the instructions provided.
- Apply a disposable probe cover to the clean, dry thermometer.
- Determine whether the measurement site will be oral, rectal or underarm.

Switching Dual Scale

NOTE: To change between Celsius and Fahrenheit readings, hold the POWER button down for 2 seconds until the "Lo" with an "M" displays and changes between °F and °C. Press the POWER button again to turn the thermometer off.

Once the desired temperature scale is selected, press the POWER button (the unit will beep) to activate the thermometer. This is the LCD display test, which ensures all parts of the display are functioning properly. The FeverVue frown face icon will also light up.

The thermometer will display the last reading taken, followed by a self-test temperature 98.6°F(37.0°C) for about 1 second. The display will show "Lo" with a flashing °F(°C) indicating the thermometer is ready for measurement.

NOTE: If the temperature is less than 89.6°F(32.00°C), the LCD displays "Lo" °F(°C). If the temperature is greater than 110°F(43.00°), the LCD displays "Hi" °F(°C).

Taking an Oral Temperature


Approximate measurement time: 10-20 seconds

- Using the diagram on the right, place the probe tip under the tongue as near as possible to the heat pocket. The diagram also shows the temperature variables by not placing the tip in the heat pocket area.

NOTE: Drinking hot or cold fluids, exercising, smoking or performing other activities may raise or lower your temperature. Therefore, it is important to relax for approximately 5 minutes with your mouth closed prior to taking a reading.

- The current temperature is shown and the °F(°C) will continue to flash during the reading.
- When the peak temperature has been reached, the °F(°C) will stop flashing. The thermometer will sound ten long beeps if the temperature is below 100°F(37.8°C). If the temperature is above 100°F(37.8°C) the unit will sound 30 short beeps and the FeverVue frown face icon will light up.
- To turn the unit off, press the ON/OFF button. The unit will automatically shut off after approximately three minutes.

NOTE: The mouth must remain closed during the measurement for the most accurate reading. Opening the mouth could result in an extended measurement time and may affect the reading.

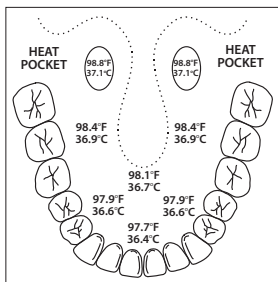Taking an Axillary Temperature (under the arm)

Approximate measurement time: 20-30 seconds

- Wipe the underarm with a dry towel.
- Place the probe tip under the arm so the tip is touching the skin with the thermometer perpendicular to the body.
- Position the arm across the chest so the probe tip is well covered by the arm. This also ensures the probe is not affected by the room's air.
- When the peak temperature has been reached, the °F(°C) will stop flashing. The thermometer will sound ten long beeps if the temperature is below 100°F(37.8°C). If the temperature is above 100°F(37.8°C) the unit will sound thirty short beeps.
- To turn the unit off, press the ON/OFF button. The unit will automatically shut off after approximately three minutes.

NOTE: Axillary temperatures are generally 1.0°F(0.5°C) lower than oral readings.

Taking a Rectal Temperature

Approximate measurement time: 10 seconds

NOTE: Rectal temperatures are an appropriate and reliable method for infants and small children. Otherwise, this method should only be used when it is impossible or impractical to take an oral or axillary temperature.

- After applying the disposable probe cover, apply a water-soluble lubricant to the tip; do not use petroleum jelly.
- Gently insert the probe a MAXIMUM of 1/2" into the rectum.
- When the peak temperature has been reached, the °F(°C) will stop flashing. The thermometer will sound ten long beeps if the temperature is below 100°F(37.8°C). If the temperature is above 100°F(37.8°C) the unit will sound thirty short beeps.
- To turn the unit off, press the ON/OFF button. The unit will automatically power down after approximately three minutes.

NOTE: Rectal temperatures are generally 1.0°F(0.5°C) higher than oral readings.

Using Disposable Probe Covers

It is recommended to use the enclosed disposable probe covers to prevent cross contamination. Additional probe covers can be purchased from your retailer or with the order form provided.

General Care

- The thermometer should only be used under the supervision of an adult.
- Do not walk, run or talk during temperature taking.
- Clean the thermometer before and after each use.
- Store the unit in the protective case when not in use.
- Do not bite the probe.
- Do not store the unit where it will be exposed to direct sunlight, dust or humidity. Avoid extreme temperatures.
- Dropping or subjecting your thermometer to strong shocks should be avoided.
- Do not attempt to disassemble the unit, except to replace the battery.

NOTE: Performance of the device may be degraded if: operated outside of stated temperature and humidity range; stored outside of stated temperature and humidity range; thermometer undergoes mechanical shock (drop); patient temperature is below ambient temperature.


Changing the Battery

The battery in your digital thermometer needs to be replaced when the "■" symbol appears in the lower right corner of the display. Replace the battery with a AG10(LR1130/LR54) alkaline or equivalent.

To replace the battery,

- Remove the battery cover from the top of the thermometer.
- Use a pointed tool to remove the battery.
- Insert the new battery with positive (+) side facing up.
- Slide the battery cover back in to place.

NOTE: Please properly dispose of the battery away from small children and heat.

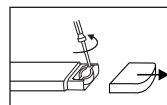Specifications

Item 15-905-000

DISPLAY:	Liquid Crystal Display
DISPLAY RESOLUTION:	0.1°F (0.1°C)
MEASUREMENT RANGE:	90.0°F (32.0°C) - 109.9°F (42.9°C)
ACCURACY:	±0.2°F; 98.0°F - 102.0°F (±0.1°C; 37.0°C - 39.0°C) ±0.3°F; 96.4°F - 98.0°F (±0.2°C; 35.8°C - 37.0°C) ±0.3°F; 102.0°F - 106.0°F (±0.2°C; 39.0°C - 41.0°C)
AVERAGE MEASUREMENT TIME:	Oral 10-20 seconds. Rectal 10 seconds. Axillary 20-30 seconds
TONE:	Thermometer beeps when activated and deactivated. After the measurement is complete, ten long beeps sound for temperatures under 100°F (37.8°C), 30 short beeps sound for temperatures above 100°F (37.8°C).
MEMORY:	Last 10 readings
BATTERY:	One, LR1130 or equivalent
BATTERY LIFE:	Approximately 1 year
AUTOMATIC SHUT-OFF:	After approximately 3 (+/- 1) minutes
AMBIENT OPERATING TEMPERATURE:	50.0°F to 104.0°F (10.0°C - 40.0°C)
STORAGE TEMPERATURE:	14.0°F - 140.0°F (-10.0°C - 60.0°C)
DIMENSIONS:	5-1/4" x 3/4" x 1/2"
WEIGHT:	0.5 ounces (13 grams)

Specifications are subject to change without notice.

LIMITED ONE-YEAR WARRANTY

The warrantor guarantees that its Thermometer will be free from manufacturing defects under normal use for one year from the original purchase date. This warranty covers only normal use, and does not apply to use in any clinical or commercial applications. This warranty does not cover batteries that may be provided with or used with the Thermometer. This warranty is voided if the Thermometer product is misused or abused in any manner.

If the Thermometer fails to operate during the first year from the time originally purchased, please call our Customer Care Help Line at 800-622-4714 and we will try to resolve the issue as soon as possible. If the warrantor determines the unit failed to operate due to a manufacturing defect, the unit will be replaced at the option of the warrantor. Replacement is the sole remedy under this limited warranty. This warranty gives you specific legal rights, which vary from state to state. As a condition of this warranty, the enclosed warranty registration card must be completed and sent to us within 10 days of the purchase date or can be filled out online at www.mabisdmi.com.

This warranty constitutes the warrantor's only responsibility and obligation to replace materials or components. We make no other express or implied warranties, arising by operation of law or otherwise, or any warranty of merchantability or fitness for a particular use or purpose whether or not the use or purpose has been disclosed to the warrantor in specifications, drawings or otherwise, and whether or not the warrantor's products are specifically designed and/or manufactured by the warrantor for the buyer's use or purposes, except for the limited warranty stated above. The warrantor will not be responsible for any indirect, incidental, special, consequential, or punitive damages or other loss, including, but not limited to, damage to or loss of other property or equipment and personal injuries, whether to purchaser or others. The warrantor shall in no event be liable to the purchaser for any amount in excess of the cost of repair and/or replacement of the unit.

CUT OUT AND MAIL TO:
WARRANTY REGISTRATION DEPARTMENT
1931 NORMAN DRIVE SOUTH
WAUKEGAN, IL 60085

OR REGISTER ONLINE AT WWW.MABISDMI.COM

Thermometer WARRANTY REGISTRATION
(MUST BE COMPLETED WITHIN 10 DAYS OF PURCHASE)

NAME: _____

ADDRESS: _____

CITY: _____ STATE: _____ ZIP: _____

PURCHASED FROM: _____

CITY: _____ STATE: _____ ZIP: _____

DATE PURCHASED: _____

MODEL PURCHASED: 15-905-000

USER INFORMATION (OPTIONAL)

MALE FEMALE AGE: _____

AGE OF YOUNGEST FAMILY MEMBER

UNDER 6 MONTHS 6 MONTHS - 1 YEAR
 1 - 3 YEARS 3 - 5 YEARS OLDER

PROBE COVER ORDER FORM

CUT OUT AND MAIL TO: 1931 NORMAN DRIVE SOUTH
WAUKEGAN, IL 60085

Please send _____ packages of disposable probe covers #15-618-000 (100 probes per package). Enclosed is \$5.00 per package plus \$3.00 for shipping and handling for each order.

PAYMENT METHOD (do not send cash): Check Money Order

Please ship to:

NAME: _____

ADDRESS: _____

CITY: _____ STATE: _____ ZIP: _____

PHONE: (____) _____

Probe covers can be purchased online at mabisdmi.com.

NOTE: This order form is for shipments within the United States and Canada.

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